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Texto revisado: 10/05/2024

INSTRUCTIONS FOR USE OF SOFT CONTACT LENSES

CE 0318 (According to Regulation 2017/745 EU)	CE 0318
Sterilized using water steam	STERILE
Medical Device	MD
Do not use if the container is damaged	
GREEN POINT (Ecoembes)	
MANUFACTURER	
BATCH	LOT
MANUFACTURING DATE	
DATE OF EXPIRY	

 iLents Care, SL
Pintores, 5
28923 Alcorcón
MADRID ESPAÑA
www.tiedra.net


HOW TO USE: Follow the steps below to put in and take out your soft contact lenses.

TO PUT IN CONTACT LENSES

- ① Before handling the lenses, wash and rinse your hands well.
- ② Place your contact lens on your fingertip, making sure it is in the correct position. Try to always proceed in the same order.
- ③ Separate and hold both eyelids with your fingers and place the lens on the cornea.
- ④ Press the contact lens lightly on the eye and slowly release the eyelids. Repeat with the other lens.

TO REMOVE CONTACT LENSES

- ⑤ To remove the lens, grasp the lower eyelid, gently move the lens downward while looking up.
- ⑥ Pinch the lens gently with the pad of your thumb and forefinger, never your fingernails, and remove

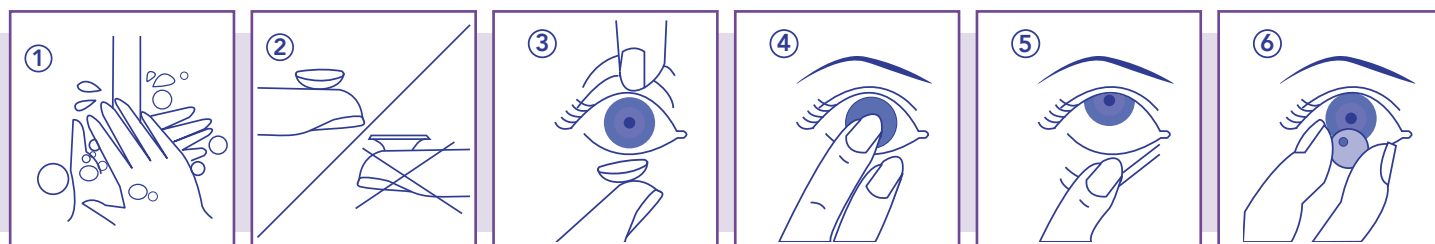INSTRUCTIONS FOR USE OF SOFT CONTACT LENSES

EN

TIEDRA

SOFT CONTACT LENSES **TIEDRA**

- These contact lenses have been recommended by your eye care professional, it is essential that you follow their recommendations at all times.
- He will tell you the most appropriate replacement time and the number of hours of use.
- Visit your eye care professional regularly (at least every 6 months) and follow their instructions.
- Never use contact lenses that have not been recommended by a professional.
- If you have any questions, consult your eye care professional.
- If these contact lenses are going to be used by a minor, the mother and/or father must pay special attention to everything indicated in these instructions for use to ensure the safe use of this contact lens by their daughter or son.
- Always follow what is indicated in these instructions for use

WARNINGS: Problems with contact lenses and lens care products may cause infections or ulcers on the cornea, leading to loss of vision. If you do not dispose of the lens holder solution after each use or if you leave the lenses

in tap, mineral or distilled water or any non-sterile solution, they may become contaminated and cause eye damage and possible loss of vision. Please review the **Important Safety Instructions** section below for the safe

use of your contact lenses.

Important security instructions. For a safe use of your contact lenses, follow the instructions below for the

Maintenance of these contact lenses:

- These soft contact lenses must be kept in a clean lens holder completely immersed in a maintenance solution that will be renewed every day.
- Always use the cleaning and care solutions that your eye care professional recommends for your contact lenses. These contact lenses can be maintained with single contact lens solution or peroxide contact lens disinfection systems. In turn, your eye care professional will be able to recommend other contact lens cleaning products according to your specific needs.
- Wash and dry your hands thoroughly before handling lenses
- **Do not use tap water, mineral water or saliva to clean the lenses or to fill the lens holder .**
- Always discard any remaining solution from the lens holder after each use.
- Do not add maintenance solution to lens holder with used solution.
- Clean, rinse and air dry the lens holder after each use.
- Place the lens holder upside down to drain excess solution as it dries.
- Replace the lens holder with a new one frequently, at least when starting a maintenance solution bottle
- Replace these contact lenses and with new ones when appropriate according to the replacement period indicated by your eye care professional.

COMPOSITION :

Trade name of the Lens you use	Composition
EXACTA 38%, EXACTA 38% Toric, EXACTA 38% MF, EXACTA 38% MF Toric, EXACTA 6, EXACTA 6T	Polymacon (62%) and Isotonic Buffered Saline (38%).
EXACTA 50%, EXACTA 50% Toric, EXACTA 50% MF, EXACTA 50% MF Toric, ESENCIA, ESENCIA Toric, PROFAQ, PROHEAL, EXACTA SK, EXACTA PL, EXACTA SK Toric, EXACTA PL Toric	Acofilcon B (51%) and Isotonic buffered saline solution (49%).
EXACTA O2, EXACTA O2 Toric, EXACTA O2 MF, EXACTA O2 MF Toric, ESENCIA O2, ESENCIA O2 Toric, PROFAQ O2, PROHEAL O2	Filcon 5B (50) (38%) and Isotonic Buffered Saline (62%)
EXACTA 59%, EXACTA 59% Toric	Acofilcon A (42%) and Isotonic buffered saline solution (58%).

HOW DO THESE CONTACT LENSES WORK?:

Trade name of the Lens you use	Explanation of operation
EXACTA 38%, EXACTA 38% Toric, EXACTA 38% MF, EXACTA 38% MF Toric, EXACTA 6T, EXACTA 59%, EXACTA 59% Toric, EXACTA 50%, EXACTA 50% Toric, EXACTA 50% MF, EXACTA 50% MF Toric, EXACTA O2, EXACTA O2 Toric, EXACTA O2 MF, EXACTA O2 MF Toric	These contact lenses are made with a geometry that is defined by the anatomical measurements of your ocular surface and the physiological requirements of your eye. With this information we generate a geometric shape in the contact lens capable of modifying the direction of the light rays, so that they focus on your retina, and thus your vision is clear.
ESENCIA, ESENCIA O2, ESENCIA Toric, ESENCIA O2 Toric	The ESENCIA® and ESENCIA® O2 lenses, in addition to the above, incorporate curves in their peripheral geometry that ensure that peripheral light rays also focus in front of your peripheral retina, which reduces the stimuli that make myopia grow more quickly.

WHAT BENEFITS DO THESE CONTACT LENSES PROVIDE?:

Trade name of the Lens you use	Benefits you get when using these Contact Lenses
EXACTA 38%, EXACTA 38% Toric, EXACTA 38% MF, EXACTA 38% MF Toric, EXACTA 6T, EXACTA 59%, EXACTA 59% Toric, EXACTA 50%, EXACTA 50% Toric, EXACTA 50% MF, EXACTA 50% MF Toric, EXACTA O2, EXACTA O2 Toric, EXACTA O2 MF, EXACTA O2 MF toric	By using these contact lenses, your vision ametropia will be compensated, so that you will achieve visual acuity (clear vision) that allows you to carry out your daily activities normally.
ESENCIA, ESENCIA O2, ESENCIA Toric, ESENCIA O2 Toric	By using these lenses, your vision ametropia will be compensated, so that you will achieve a visual acuity (clear vision) that allows you to carry out your daily activities normally. In turn, the growth rate of the progression of myopia will be slowed down, which will reduce the final graduation of myopia.

INTENDED PURPOSE :

Trade name of the lens you use	Intended Purpose
EXACTA 38%, EXACTA 38% Toric, EXACTA 38% MF, EXACTA 38% MF Toric, EXACTA 6T, EXACTA 59%, EXACTA 59% Toric, EXACTA 50%, EXACTA 50% Toric, EXACTA 50% MF, EXACTA 50% MF Toric, EXACTA O2, EXACTA O2 Toric, EXACTA O2 MF, EXACTA O2 MF Toric	Compensation for ametropia of vision.
ESENCIA, ESENCIA O2, ESENCIA Toric, ESENCIA O2 Toric	Controlling the progression of myopia.

INDICATIONS (USES):

Trade name of the Lens you use	Intended Purpose
EXACTA 38% EXACTA 50% EXACTA 59% EXACTA O2 EXACTA 6T	Myopia and hyperopia compensation in adults and children.
EXACTA 38% Toric EXACTA 50% Toric EXACTA 59% Toric EXACTA O2 Toric	Myopia, hyperopia and astigmatism compensation in adults and children.
EXACTA 38% MF EXACTA 50% MF EXACTA O2 MF	Presbyopia compensation for adults.
EXACTA 38% MF Toric EXACTA 50% MF Toric EXACTA O2 MF Toric	Presbyopia and astigmatism compensation in adults.
ESENCIA ESENCIA O2	Compensation and Control of myopia progression in paediatric and adolescent users.
ESENCIA Toric ESENCIA O2 Toric	Myopia and astigmatism compensation. Control of myopia progression in paediatric and adolescent wearers.

CONTRAINDICATIONS (REASONS NOT TO USE THIS PRODUCT): You should not wear contact lenses in the following situations. If you have any questions, consult your optician-optometrist or ophthalmologist:

1. When there are ocular pathologies, such as: ocular infections, conjunctivitis or inflammation of the conjunctiva, keratitis or corneal inflammation, dry eye, acute inflammation of the anterior segment of the eye (iritis), reduced corneal sensitivity, ptosis, chronic ocular allergy .
2. Systemic pathologies that may affect the corneal surface.
3. Some medications, especially for the eyes, are contraindicated in the use of contact lenses. If this is your situation, consult your optician-optometrist or ophthalmologist.
4. Users who, for whatever reason, are not able to maintain or handle contact lenses properly.

PRECAUTIONS AND THINGS TO KEEP IN MIND WHEN USING CONTACT LENSES:

- Do not use the lenses if the original vial is broken, deteriorated or opened. In this case, go to the center where the contact lenses have been fitted.
- The original outer container should be stored until the last lens in its contents has been used.
- If you do not wear contact lenses for a few days, ask your optician-optometrist for advice on their correct storage.
- Do not practice water sports with the lenses on without wearing goggles that do not let the water through.
- Put your contact lenses in before putting on eye makeup and remove them before removing makeup. Use cosmetic products recommended for contact lens wearers. Makeups can color your lenses.
- Avoid contact of the lenses with soaps, creams, shampoos, lacquers, or other cosmetics as they can irritate your eyes. Close your eyes if you use hairspray, as it can damage your glasses.
- Soft contact lenses are for daytime wear. Remember to remove contact lenses before bed. On some occasions, the ophthalmologist can prescribe the continued use of soft contact lenses, in this particular case, follow what is indicated by your Ophthalmologist.
- Never use totally or partially dehydrated lenses.
- Keep contact lenses and maintenance products out of the reach of children.
- You should not use the product after the expiration date stated on the carton and bottle.

ADVERSE REACTIONS (PROBLEMS AND WHAT TO DO): During the use of contact lenses and their solutions for maintenance, you may experience the following reactions: stinging, burning or itching (irritation), reduced comfort compared to the first time you put the lenses on, sensation of something in the eye (a foreign body, a scratch in a certain area), excessive tearing of the eye, unusual discharge from the eye, red eyes, corneal oedema, corneal neovascularization, abnormalities of the tarsal conjunctiva, iritis, conjunctivitis, decreased visual acuity (reduced vision), blurred vision, flashes or halos around objects, sensitivity to light (photophobia) or dry eyes.

If you experience any of these reactions: **Remove lenses immediately and**

- If the discomfort or problem goes away, carefully examine the lenses closely.
- If the lens is damaged, do not put it back in your eye. Put the lens in its case and contact your eye care professional.
- If the lens is dirty, has a tab or a foreign body, or if the problem goes away and the lens is undamaged and has been thoroughly cleaned, rinsed, and disinfected, re-insert the lens.
- **If the problem persists, remove the lenses and see your eye care professional immediately**

NOTIFICATION OF ADVERSE EFFECTS:

You can notify us of medical incidents or adverse reactions at the web address:
<http://ae.tiedra.net/>.

You can also notify the Health Authorities of the country in which you are.

Disposal: Dispose of the product in the normal way. Try to recycle cardboard and plastic according to what is indicated in your city

Presentation: the product is presented in different formats. The lot number and expiration date are listed on the product bottle.

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